



SpectroElectroChem Suite

USER MANUAL

Version 5.5.3

Open-source software for Raman, SERS, absorption and fluorescence spectroelectrochemistry, and RRDE data analysis

Author	Prof. Dr. Achim Habekost
Release	v5.5.3
Concept DOI	https://doi.org/10.5281/zenodo.21283231
Repository	https://github.com/Achim-Habekost/SpectroElectroChem-Suite
License	MIT License

Parts of the source code were developed with the assistance from OpenAI ChatGPT and were subsequently reviewed, scientifically validated, modified and substantially extended by the author.

Document status and scope

This manual describes SpectroElectroChem Suite version 5.5.3. It is intended for scientists, students, and technical staff who wish to analyze their own comma-separated value (CSV) files without modifying the Python source code. The software provides graphical workflows for the processing and visualization of Raman, SERS, absorption, fluorescence and RRDE data for various types of spectra and voltammograms, e.g. for smoothing, background subtraction and creating waterfall, surface and other diagrams.

Important: SpectroElectroChem Suite is now available in two forms. (1) A ready-to-run Windows program: a self-contained executable (SpectroElectroChem_Suite.exe) that can be started directly, plus an optional Windows installer that adds Start-menu and desktop shortcuts. Neither requires a separate Python installation. (2) The Python source distribution, which is started with the included batch files and requires Python 3.10 or newer. Both forms provide the same four analysis modules (Raman, SERS/Raman voltammograms, absorp-to-/fluorovoltammograms and RRDE).

Table of contents

1. Overview
2. Obtaining the software
3. System requirements
4. Installation on Windows
5. Starting and closing the application
6. Input-file specifications
7. Raman spectrum analysis
8. Raman/SERS voltammograms
9. Absorption and fluorescence voltammograms
- 10 RRDE Analysis
- 11 Signal processing
- 12 Graphical output
- 13 Excel export
- 14 Recommended scientific workflow
- 15 Troubleshooting
- 16 Reproducibility
17. Citation and licensing
- Appendix A. CSV templates
- Appendix B. Output-file glossary

Appendix C. Quick-start checklist

Appendix D. Typical RRDE Parameters

1. Overview

1.1 Purpose

SpectroElectroChem Suite is a modular desktop application for processing, visualizing and exporting spectroscopic and spectro-electrochemical data. It is designed to reduce the amount of manual spreadsheet work required for multi-spectrum data sets while preserving transparent numerical outputs for independent inspection.

1.2 Supported analytical workflows

- Single Raman spectra or multiple Raman spectra stored in columns: smoothing, baseline correction, peak analysis and waterfall visualization.
- Raman and SERS voltammograms: potential-resolved Raman intensity matrices.
- Absorptovoltammograms: potential-resolved absorption spectra.
- Fluorovoltammograms: potential-resolved fluorescence spectra.
- RRDE voltammograms: rotation-rate-resolved disk and ring currents for the rotating ring–disk electrode.

1.3 Principal outputs

- Processed Excel workbooks containing numerical results and metadata.
- Static PNG and PDF figures for reports and publications.
- Interactive HTML figures for rotation, zooming and visual inspection.
- 3D surfaces, heatmaps, contour plots and waterfall diagrams.
- Raw, smoothed, baseline-corrected and vertically offset data, depending on the selected module.
- Levich, Koutecký–Levich, Tafel and Butler–Volmer analyses.

1.4 Intended use and limitations

The software is intended as a scientific analysis and visualization aid. It does not replace validation of instrumental data, inspection of raw spectra, calibration, uncertainty analysis or expert interpretation. Parameters such as the smoothing window, baseline model, peak prominence and waterfall offset must be selected with reference to the spectral resolution, the signal-to-noise ratio and the specific scientific question.

2. Obtaining the software

2.1 Zenodo archive

The archived version can be accessed through the persistent DOI:

<https://doi.org/10.5281/zenodo.21283231>

On the Zenodo record page, download the ZIP archive associated with version 5.5.3. Zenodo provides a permanently archived snapshot together with citation metadata.

2.2 GitHub repository

The development repository and release page are available at:

<https://github.com/Achim-Habekost/SpectroElectroChem-Suite>

Use GitHub for source inspection, issue reporting and future releases. Use the Zenodo DOI when citing the archived version used in a scientific study.

2.3 After downloading

- Save the downloaded ZIP archive to a local folder.
- Right-click the archive and choose Extract All.
- Open the extracted project folder. Do not run batch files directly from inside the ZIP archive.

Security notice: Windows may mark downloaded batch files as downloaded from the Internet. Inspect the files and source code before execution, and use only the official GitHub or Zenodo release.

3. System requirements

Component	Minimum / required	Recommended
Operating system	Windows 10 or 11	64-bit Windows 11
Python	Python 3.10 or newer	Current stable 64-bit Python
Memory	4 GB RAM	8 GB or more for large matrices
Storage	Approx. 1 GB plus output files	Several GB for large Excel/HTML exports
Internet	Needed initially for package installation	Needed to load Plotly from a CDN in some HTML outputs
Spreadsheet viewer	Optional	Microsoft Excel or LibreOffice Calc
Web browser	Modern browser	Firefox, Edge or Chrome

3.1 Required Python packages

The release requirements file lists the packages used by the different modules. The principal packages are:

- NumPy – numerical arrays and matrix operations.
- pandas – CSV input and tabular data handling.
- SciPy – Savitzky–Golay filtering, peak processing and numerical routines.
- pybaselines – baseline-estimation algorithms, where used.
- Matplotlib – static plots and PDF/PNG export.
- Plotly – interactive HTML visualizations.
- openpyxl – Excel workbook creation.
- PySide6 and/or Tkinter – graphical user interfaces, depending on the module implementation.

4. Installation on Windows

Revision (v5.5.3): The self-contained Windows program — the ready-to-run executable and the optional installer described in Section 4.5 — is the recommended installation route and requires no separate Python installation. Installation from the Python source distribution (Sections 4.1–4.4) is provided as an alternative for advanced users and developers and is summarised in Appendix E.

4.1 Install Python

Alternative (advanced) route: the following Python-source steps (Sections 4.1–4.4) are only required if you do not use the executable or installer from Section 4.5. A condensed version is given in Appendix E.

- Download Python 3.10 or newer from the official Python website.
- During installation, enable Add Python to PATH if that option is shown.
- Complete the installation and restart Windows, or at least reopen File Explorer and Command Prompt.

4.2 Install required packages

In the extracted program folder, double-click:

```
Install_required_Python_packages.bat
```

A console window will open and install the packages listed in requirements.txt. Internet access is normally required during this step. The step is generally needed only once for a given Python environment.

4.3 Test the installation

Run:

```
Test_Installation.bat
```

The test verifies that Python and the principal dependencies can be imported. If the console window closes immediately, launch the batch file from Command Prompt to retain the error message.

4.4 Manual installation alternative

Alternatively, open Command Prompt in the program folder and run:

```
py -m pip install -r requirements.txt
```

If the py launcher is unavailable, try:

```
python -m pip install -r requirements.txt
```

Multiple Python installations: If package installation succeeds but the application still reports missing modules, the batch files may be using a different Python installation from the one in which the packages were installed. Use “py -0p” in Command Prompt to list the installed Python interpreters.

4.5 Installation using the Windows executable or installer

From version 5.5.3, the suite is also distributed as a ready-to-run Windows program that does not require a separate Python installation. Two variants are provided.

Portable version. Extract the distributed program folder and double-click the executable:

```
SpectroElectroChem_Suite.exe
```

All required libraries are bundled inside the folder. Keep the executable together with the accompanying files and sub-folders; do not move the .exe out of its folder.

Windows installer (setup). Alternatively, run the setup program (for example SpectroElectroChem_Suite_Setup.exe). The installer copies the program into the standard applications folder and creates Start-menu and desktop shortcuts. The application can then be started like any other installed Windows program.

The executable and the installer are built from the same source code as the Python distribution and provide the identical four analysis modules. No Python interpreter, batch file or package installation is needed for these variants.

Security notice: Because the executable is newly published and not code-signed, Windows SmartScreen or the antivirus software may display a warning on first start. Verify that the file was obtained from the official GitHub or Zenodo release, then choose “More info” and “Run anyway” if you trust the source.

5. Starting and closing the application

5.1 Standard start

Double-click:

```
Start_SpectroElectroChem_Suite.bat
```

The launcher adds the src folder to the Python search path and opens the central application window. Depending on the release build, the main window offers module buttons or menu entries

for Raman spectrum analysis, Raman/SERS voltammograms, and absorption/fluorescence voltammograms.

5.2 Start from Command Prompt

```
cd "C:\path\to\SpectroElectroChem_Suite"
python main.py
```

5.3 Starting from the Windows executable

If the ready-to-run program or the installed version is used, start the suite directly:

- Portable version: double-click SpectroElectroChem_Suite.exe inside the program folder.
- Installed version: use the Start-menu entry or the desktop shortcut created by the installer.

The central application window opens with buttons or menu entries for the four modules. No Command Prompt, batch file or Python installation is required in this case.

5.4 Closing

Close module windows using their normal Close or X button. Allow ongoing Excel or HTML export operations to finish before closing. Interactive HTML plots remain available as independent files after the application has been closed.

6. Input-file specifications

6.1 General CSV rules

- Use a true CSV text file, not an .xlsx file merely renamed to .csv.
- Decimal points and decimal commas are accepted by many import paths, but a decimal point is recommended for portability.
- Do not mix text comments into the numerical matrix unless the module explicitly supports metadata rows.
- Avoid merged cells, hidden rows and formulas; export numerical values only.
- Remove completely empty trailing rows and columns where practical.
- Keep a copy of the unmodified raw instrument export.

6.2 Matrix format for Raman/SERS voltammograms

The first row contains the potentials, the first column contains the Raman shifts (wavenumbers), and all remaining cells contain Raman intensities.

	-0.40	-0.30	-0.20	-0.10
300	120	125	131	140
305	118	124	130	139
310	119	126	134	143

	-0.40	-0.30	-0.20	-0.10
315	123	129	138	148

6.3 Matrix format for absorption and fluorescence voltammograms

The first row contains the potentials and the first column contains the wavelengths. The remaining cells contain absorbance, fluorescence intensity or another measured optical response. If there is only one absorbance or fluorescence spectrum, the potential row is missing.

	0.00	0.10	0.20	0.30
400	0.021	0.028	0.039	0.051
405	0.023	0.031	0.043	0.056
410	0.026	0.035	0.048	0.062
415	0.030	0.040	0.054	0.069

6.4 Raman spectrum table format

For Raman spectrum analysis, the first numerical column contains the Raman shift and the remaining columns contain one or more intensity spectra. Header labels may be used to identify individual spectra at various potentials. The exact import options are displayed by the module.

6.5 Axis order

Some modules sort the potential and spectral axes automatically. Nevertheless, ascending, unique numerical axis values are recommended. Duplicate potentials or duplicate wavelengths/wavenumbers should be resolved before analysis unless they represent intentional replicate measurements.

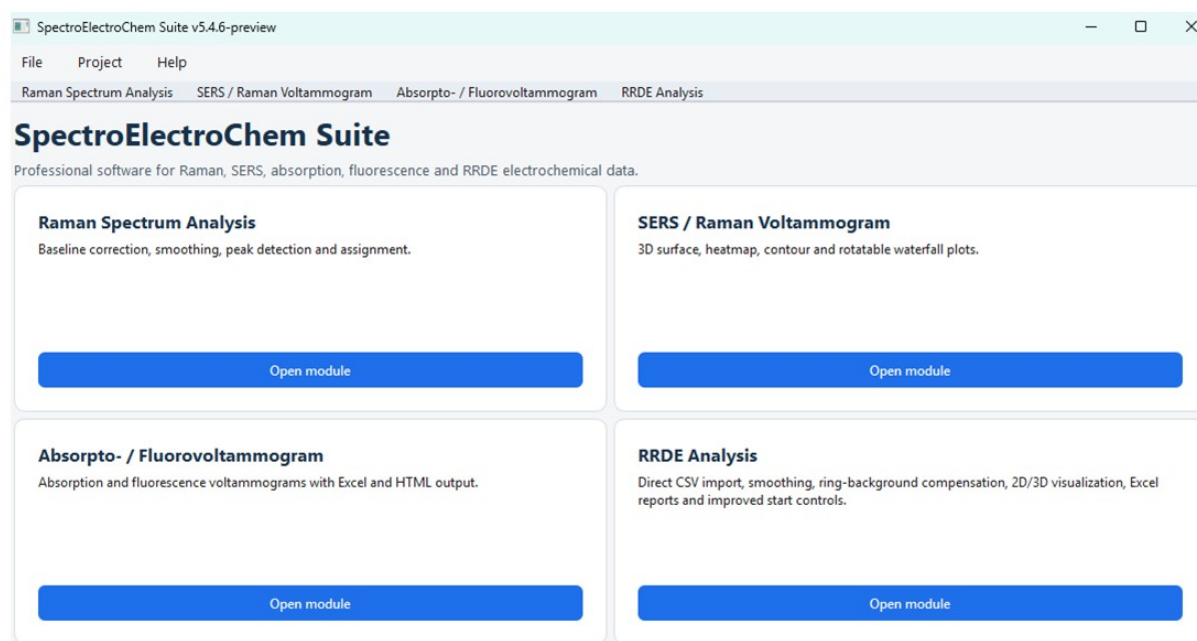


Figure 1: Initial input

7. Raman spectrum analysis

7.1 Purpose

This module processes individual Raman spectra or collections of spectra stored in columns. It can generate raw, smoothed, baseline-corrected and combined data sets, locate peaks, calculate peak parameters and create waterfall plots.

7.2 Typical workflow

- Open the Raman Spectrum Analysis module.
- Select the CSV input file.
- Enter or confirm the Raman shift range.
- Choose the baseline correction parameters.
- Choose the Savitzky–Golay smoothing parameters.
- Set the peak-finding parameters, if peak analysis is required.
- Select an optional reference peak list, if identification is required.
- Set the maximum Raman intensity for plotting, if a fixed scale is desired.
- Set the waterfall vertical offset and plot orientation, where available.
- Run the analysis and inspect the generated control plots before accepting the results.

7.3 Baseline correction

Baseline correction estimates the slowly varying background and subtracts it from the measured spectrum. A visually smooth baseline-corrected spectrum is not sufficient evidence that the baseline

is physically appropriate. Depending on the current plugin, the module may provide polynomial or penalized-baseline methods. Baseline correction should remove broad fluorescence or instrumental background (e.g. scattering) without distorting genuine broad Raman bands.

Scientific check: Always inspect the baseline overlay. A visually smooth baseline-corrected spectrum candidate assignment is not sufficient evidence that the baseline is physically appropriate. Check low-intensity regions, broad bands and the spectrum edges.

7.4 Savitzky–Golay smoothing

Savitzky–Golay smoothing performs a local polynomial fit. The window length must be odd and larger than the polynomial order. A window that is too large can reduce peak heights, merge neighbouring bands or alter the FWHM. Select a window consistent with the sampling interval and the expected peak width.

7.5 Peak analysis

Peak finding typically uses prominence and minimum-distance criteria. The output may include peak position, intensity and full width at half maximum (FWHM). Peak identification against a reference list is tolerance-based and must be treated as a tentative assignment, not as definitive chemical proof.

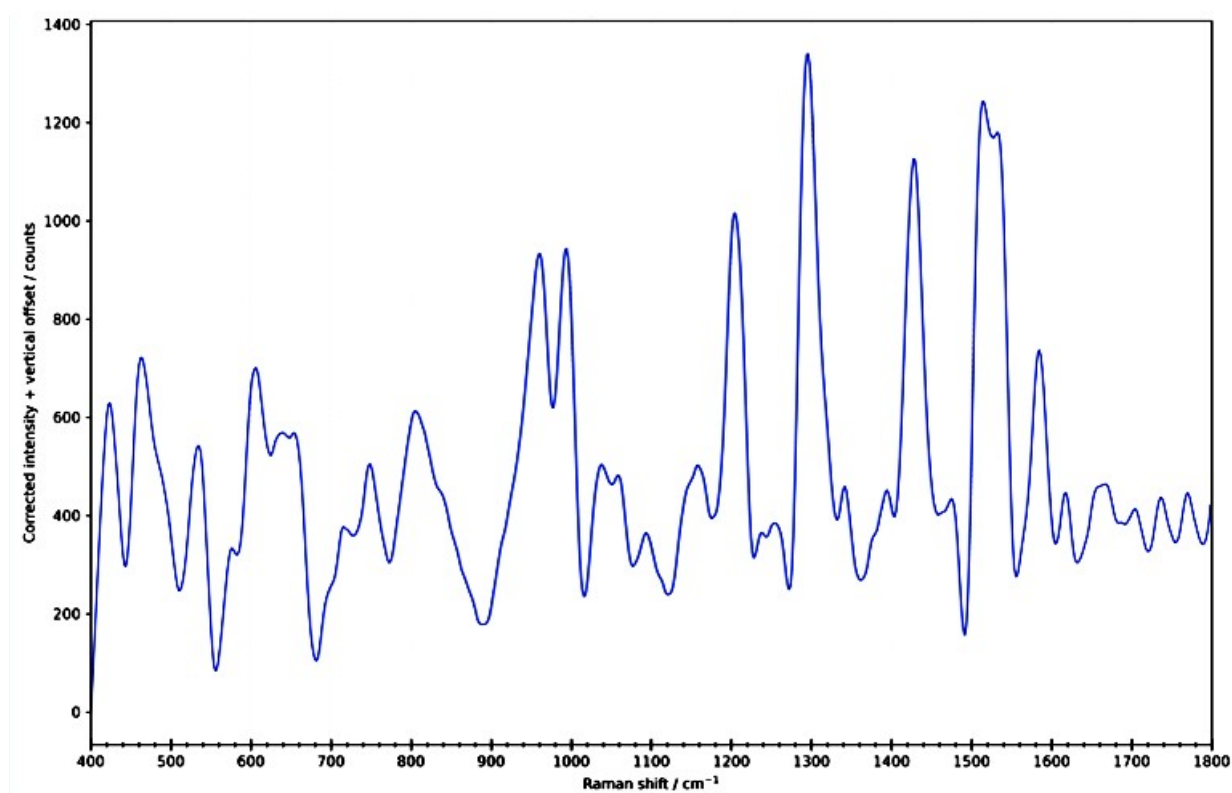


Figure 2: Smoothed and baseline-corrected Raman spectrum (rubrene).

7.6 Raman Excel workbook

Depending on the software version the Excel workbook typically contains the following worksheets:

- Original spectra
- Smoothed-only spectra
- Baseline-corrected-only spectra
- Baseline-corrected and smoothed spectra
- Calculated baselines
- 2D waterfall data
- Peak parameters
- Peak identification results
- Reference data used
- Summary/interpretation

8. Raman/SERS voltammograms

8.1 Purpose

The Raman/SERS voltammogram module processes a spectral intensity matrix as a function of Raman shift and electrode potential. It is suitable for potential-dependent Raman or surface-enhanced Raman data.

8.2 Loading and range selection

- Select the CSV matrix file.
- Select an output folder, or leave the output-folder field empty to use the input-file folder.
- Use the range-reading function, if available, to display the detected Raman-shift and potential ranges.
- Enter the desired lower and upper Raman-shift limits.

8.3 Smoothing options

- None – preserves the measured matrix.
- Moving average – simple smoothing over a selected number of points.
- Savitzky–Golay – polynomial smoothing with a selected odd window length and polynomial order.

8.4 Intensity scaling

- Raw / linear – preserves the original plotted magnitude.
- Clip at the 99th percentile + log1p reduces the influence of extreme spikes on the color scale.

- log1p – compresses the dynamic range after shifting negative values, if necessary.
- Clip99 + log1p – combines spike suppression and logarithmic compression.

Interpretation: Scaling modes modify the visualization matrix and may also be exported as plotted values. They do not convert arbitrary intensity units into quantitative concentrations.

8.5 Baseline correction and peak parameters

- SERS spectra usually require little or no baseline correction.

Peak No	Position (cm ⁻¹)	Intensity	FWHM (cm ⁻¹)
1	443.572	3237.016	70.73005
2	485.158	2005.668	13.27788
3	545.501	104.1584	3.79017
4	607.772	423.8645	59.87418
5	642.605	179.7476	6.285982
6	708.867	-260.501	2.38847
7	763.801	989.6835	41.52643
8	807.793	225.7842	3.46681
9	853.913	-26.1113	9.527231
10	879.333	-189.92	3.532601
11	909.649	27.43816	5.6194
12	954.741	271.4507	79.89512
13	991.971	247.2508	15.23369
14	1067.933	55.90622	56.92809
15	1171.156	1220.993	47.16028
16	1220.683	11.66072	9.632227
17	1246.398	-48.9973	6.112788
18	1299.668	530.7832	35.58292
19	1331.767	-33.6803	3.2529
20	1395.234	2346.762	33.20733
21	1457.745	178.5831	3.400773
22	1486.447	410.5359	60.54273
23	1547.589	-144.147	35.24608
24	1584.265	-229.842	11.93203
25	1618.48	2967.378	34.49392
26	1671.359	157.6585	34.36061

8.6 Manual intensity limits

When automatic intensity scaling is disabled, enter minimum and maximum intensity values. Interactive 3D figures may also include buttons for automatic Z scaling, manual Z scaling and percentile-based ranges.

8.7 Waterfall and surface plots

The waterfall plot displays each spectrum with a vertical displacement. Use the vertical offset to prevent overlap while preserving relative peak patterns. Interactive plots can be rotated; static plots are intended for publications and reports. Avoid offsets so large that potential-dependent trends become visually exaggerated.

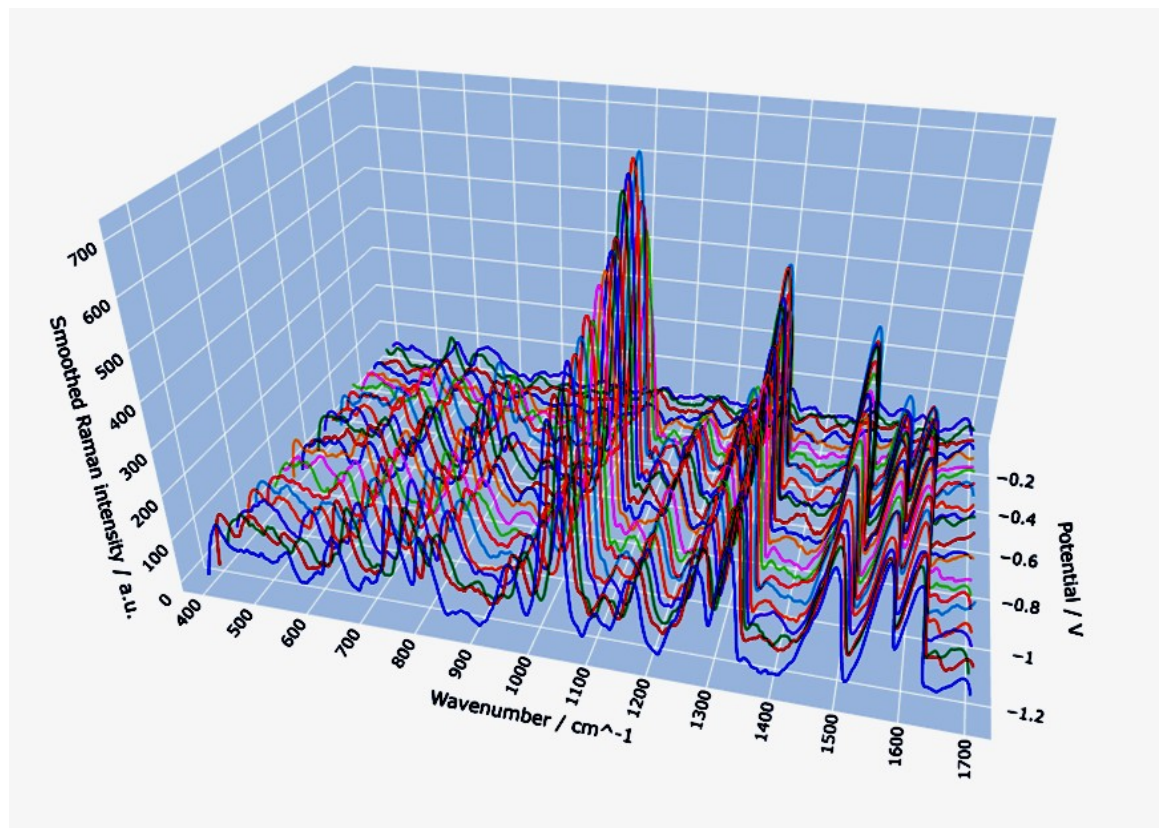


Figure 3: Interactive 3D waterfall plot of methylene blue Raman spectra recorded at different electrode potentials.

8.8 Expected outputs

- Processed Excel workbook
- Interactive 3D surface (HTML)
- Interactive heatmap (HTML)
- Interactive contour plot (HTML)
- Interactive rotatable waterfall plot (HTML)
- Optional raw-versus-smoothed comparison waterfall
- Static PNG/PDF plots, where enabled

8.9 Excel sheets

The exported workbook typically contains Metadata, Raw_Intensity_Matrix, smoothed_Intensity_Matrix, Plotted_Intensity_Matrix and the corresponding long-format tables.

Current waterfall-enabled releases may also include shifted values, unshifted values and offset tables for custom plotting in Excel.

9. Absorpto- and fluorovoltammograms

9.1 Purpose

This combined module analyzes potential-resolved absorption or fluorescence matrices. The data structure is analogous to that of the Raman voltammogram module, but the spectral axis is wavelength rather than Raman shift, and the response values are absorbance or fluorescence intensity.

9.2 Select the analysis mode

Choose Absorption or Fluorescence so that titles, axis labels and output metadata match the selected measurement type. The underlying matrix orientation remains the same: potentials in the first row, wavelengths in the first column.

9.3 Wavelength range

Select only the spectral interval needed for interpretation. Restricting the range can reduce file size and improve visual resolution. Ensure that the selected interval includes the complete band shape when comparing peak positions or widths.

9.4 Waterfall data for Excel

In addition to the graphical waterfall output, the module exports numerical waterfall values for user-defined Excel plots. The shifted matrix contains the vertical offsets, the unshifted matrix preserves the processed spectral values, and a separate offset table documents the applied vertical offset.

9.5 Typical outputs

- Processed absorption or fluorescence matrix
- Interactive surface, heatmap and contour files
- Interactive and/or static waterfall plot

- Excel workbook containing raw/processed data, metadata and waterfall.

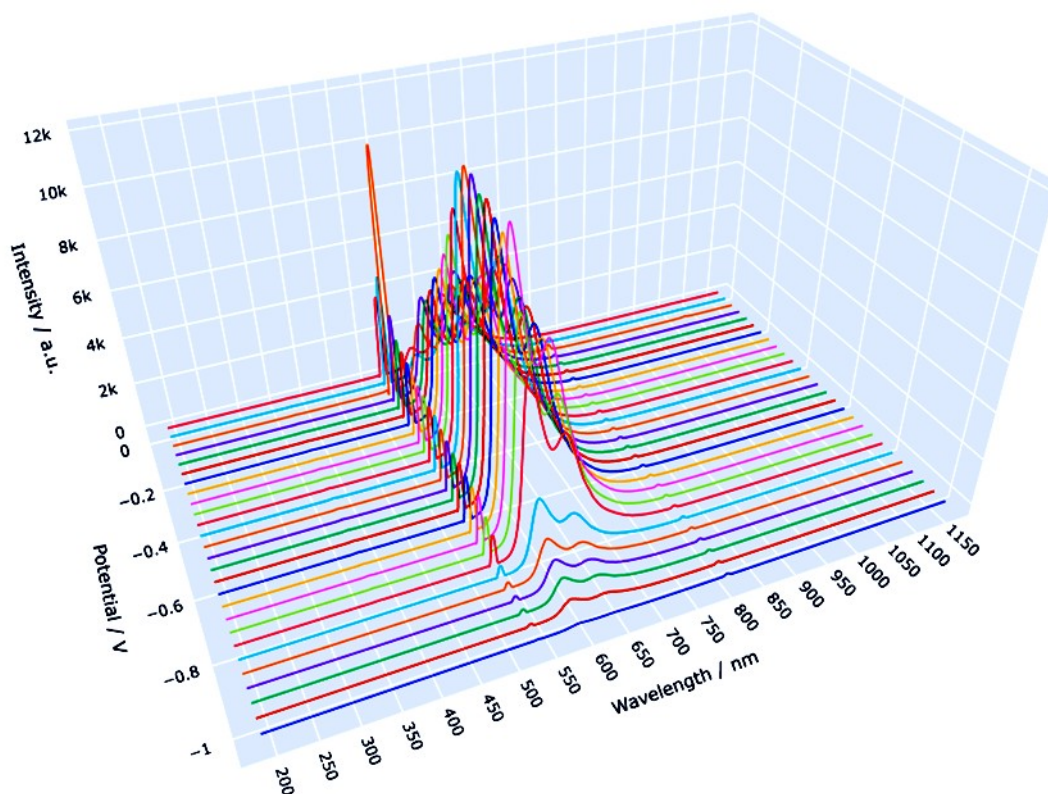


Figure 5: Fluorovoltammogram of rubrene at different potentials.

10. Rotating Ring–Disk Electrode (RRDE) Analysis

10.1 Introduction

The rotating ring–disk electrode (RRDE) is a well-established electrochemical technique for investigating electrochemical reaction mechanisms and soluble intermediates. In contrast to conventional Rotating Disk Electrode (RDE) experiments, the RRDE combines a central disk electrode with a concentric ring electrode. During an experiment, both electrodes rotate simultaneously at a precisely controlled speed.

Species generated at the disk electrode are transported by convection to the surrounding ring electrode, where they can be detected electrochemically. This enables simultaneous investigation of the primary electrochemical reaction occurring at the disk and the fate of soluble intermediates.

Typical applications include:

- Oxygen reduction reaction (ORR)
- Hydrogen peroxide detection
- Evaluation of catalyst selectivity
- Metal deposition and dissolution
- Electrochemical reaction mechanism studies

- Homogeneous redox mediators

Compared with conventional RDE measurements, RRDE experiments provide considerably more mechanistic information because both the main reaction and intermediate products can be monitored simultaneously.

10.2 Principle of the RRDE

The RRDE consists of two independently controlled working electrodes:

- **Disk electrode (WORKING ELECTRODE1):** where the primary electrochemical reaction occurs.
- **Ring electrode (WORKING ELECTRODE2):** detects products formed at the disk.

Both electrodes rotate with the same angular velocity. Rotation generates a well-defined hydrodynamic flow that transports dissolved species from the disk to the ring.

Only a fraction of the products generated at the disk reaches the ring. This fraction is called the **collection efficiency (N)**.

Typical collection efficiencies range from approximately 0.20 to 0.40, depending on the electrode geometry.

For the Metrohm Pt/GC RRDE (3.109.4080), the nominal collection efficiency is approximately **25 %**.

RRDE-Analysis
Disk- und Ringströme: Glättung, 2D, getrennte 3D-Darstellung, gemeinsame 3D-Darstellung, Excel-Report und optionale Levich/Koutecky-Levich-Analyse.

Eingabe und Ausgabe

Messdatei (CSV): Durchsuchen
Dateistatus: Noch keine CSV-Datei ausgewählt. CSV-Struktur ...
Ausgabeordner: Durchsuchen
Messungsname:

AUSWERTUNG STARTEN (oder Taste F5)

CSV file wird direkt ohne Prüfdialog ausgewertet.

Glättung

☒ Savitzky-Golay-Glättung anwenden
Fensterbreite / Punkte: Polynomgrad:

Darstellung

Disk current-Einheit: Ring current-Einheit:
☒ Rohkurven zusätzlich schwach darstellen ☐ Potentialachse umkehren
☒ Ergebnisordner nach Abschluss öffnen
Ring current-Skalierung im gemeinsamen 3D-Diagramm:
☒ Automatisch ☐ Manuell Multiplikationsfaktor:

Untergrundkompensation

☐ Untergrund kompensieren
☒ Mittelwert im Potentialbereich
Von / V: Bis / V:
☐ Manuellen konstanten Offset abziehen

Figure 6: The first part of the input for the RRDE analysis.

10.3 Data Import

SpectroElectroChem Suite imports RRDE data directly from CSV files exported by NOVA.

The expected data structure is

Column	Content
1	Potential
2	Disk current (Rotation 1)
3	Ring current (Rotation 1)
4	Disk current (Rotation 2)
5	Ring current (Rotation 2)
...	Additional rotation speeds

The first row contains the rotation speed for each measurement pair.

Example:

Potential (V)	Disk 500 rpm	Ring 500 rpm	Disk 1000 rpm	Ring 1000 rpm
0.20	...	...	...	...

After import, the software automatically identifies

- rotation speed,
- disk current,
- ring current,
- measurement order.

No manual column assignment is required, which minimizes user interaction and reduces the risk of assignment errors.

10.4 Main RRDE Window

The RRDE module provides a graphical interface for importing, processing and evaluating RRDE measurements.

The main window contains the following sections:

Data Import

- Load CSV file
- Automatic recognition of rotation speeds
- Automatic assignment of disk and ring currents

Background Correction

The user may choose between

- no correction,
- manual background correction,
- potential-interpolated subtraction of a separate N_2 /background RRDE measurement.

In version 5.5.3, background correction can be applied to both disk and ring currents. When a separate N₂/background RRDE measurement is used, the software matches the background data to the O₂ measurement by rotation rate and interpolates the background currents onto the potential grid of the corresponding O₂ measurement before subtraction.

Collection Efficiency

The collection efficiency N may either

- be entered manually,
- or taken from the electrode manufacturer's specifications.

Typical values range between 0.20 and 0.40.

10.5 Background Correction

In RRDE experiments, measured disk and ring currents may contain background contributions that are not part of the electrochemical process of interest. These contributions can influence quantitative evaluation, particularly when the ring current is small.

Examples include

- capacitive currents,
- oxidation of impurities,
- oxygen-independent reactions,
- electrode surface processes.

SpectroElectroChem Suite therefore provides two background-correction approaches.

Manual Background

A constant background current can be subtracted when the background is approximately independent of potential. This option is mainly intended for simple exploratory corrections.

For potential-dependent backgrounds, a separately measured N₂/background data set is preferable.

Nitrogen Background Measurement

A separate RRDE background measurement recorded under N₂ (or under otherwise appropriate background conditions) can be imported in addition to the O₂ measurement. The background file should contain the same rotation rates and corresponding disk/ring current pairs as the measurement to be corrected.

For each rotation rate, the software first checks the correspondence of the measurement series. Because the potential grids of the O₂ and N₂ measurements may differ slightly, the background currents are linearly interpolated onto the potential grid of the O₂ measurement before subtraction:

$$IDisk,corr(E) = IDisk,O_2(E) - IDisk,N_2,interp(E)$$

$$IRing,corr(E) = IRing,O_2(E) - IRing,N_2,interp(E)$$

The corrected currents are then used for subsequent RRDE evaluation. The program validates the point count, scan direction, potential range and rotation-rate matching and reports incompatible data sets instead of silently subtracting mismatched curves.

Matching and interpolation

Background measurements are matched to the O₂ measurement by rotation rate. For every matched pair, the N₂/background disk and ring currents are interpolated as functions of potential onto the potential values of the O₂ curve. Subtraction is therefore performed at identical potentials even when the two measurements were recorded on slightly different potential grids.

Excel documentation

The RRDE Excel report documents the complete correction procedure in a dedicated background-correction worksheet. It includes the original O₂ currents, the interpolated N₂/background currents, the background-corrected disk and ring currents and, where smoothing is enabled, the smoothed-corrected currents. This preserves the numerical steps of the correction for independent inspection and reproducibility.

Scientific check: Background subtraction is only meaningful when O₂ and N₂/background measurements were recorded under comparable experimental conditions. In particular, electrode state, electrolyte composition, scan direction, rotation rates and potential range should be consistent.

10.6 Collection Efficiency

The collection efficiency N describes the fraction of reaction products formed at the disk that reaches the ring electrode.

It depends exclusively on the RRDE geometry.

Typical values are

- Metrohm Pt/Pt RRDE: approximately 0.25
- Metrohm Pt/GC RRDE: approximately 0.25

The user-defined collection efficiency is applied in all RRDE calculations. The ratio $I_{\text{Ring}}/I_{\text{Disk}}$ may be used only under suitable experimental conditions to estimate the effective collection efficiency and should not be interpreted as its general definition. Only for calibration and subsequently applies the user-defined collection efficiency in all calculations.

An incorrect value of N directly affects

- hydrogen peroxide yield,
- electron-transfer number,
- kinetic analysis.

Therefore the manufacturer-recommended value should normally be used.

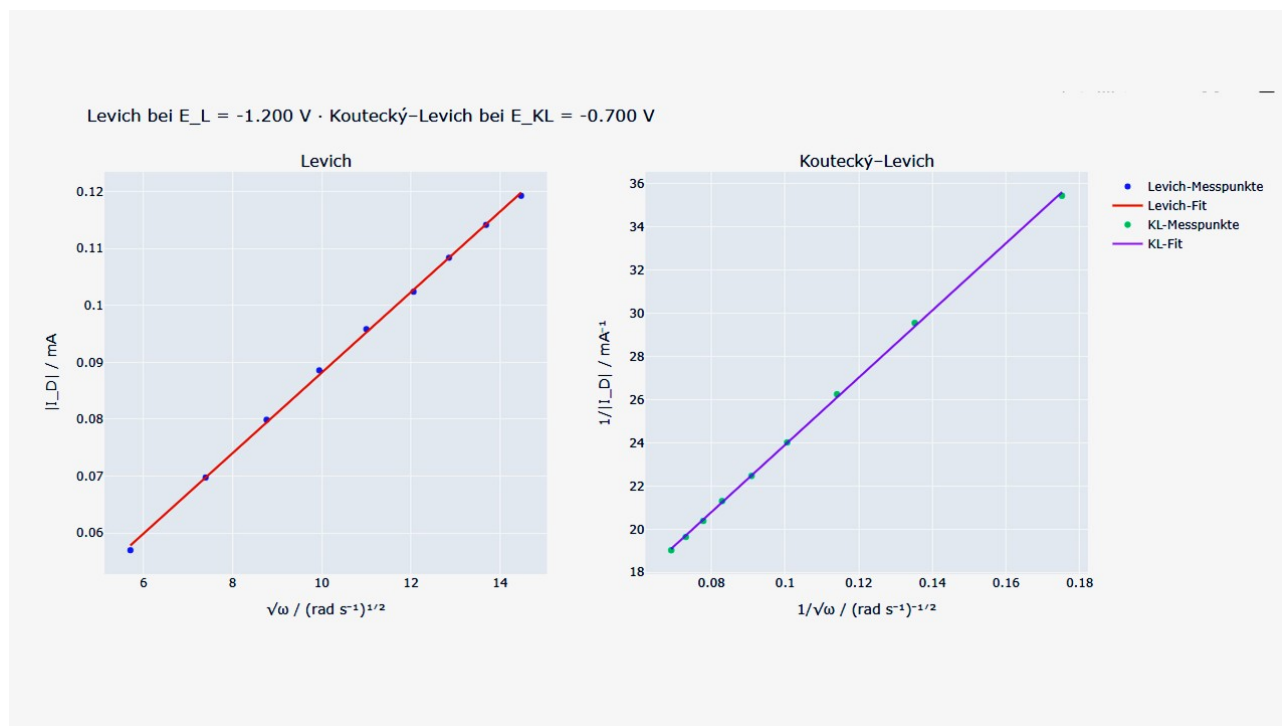


Figure 7: Levich (left) and Koutecký-Levich (right) graphs.

10.7 Electron-Transfer Number

One of the most important RRDE parameters is the apparent number of transferred electrons.

SpectroElectroChem Suite calculates

$$n = 4 \cdot I_D / (I_D + I_R/N)$$

where

- I_D = disk current
- I_R = ring current
- N = collection efficiency.

Typical interpretation:

n	Interpretation
≈ 4	Direct four-electron reduction
3–4	Mixed mechanism
≈ 2	Hydrogen peroxide formation dominates

For oxygen reduction, values close to four indicate efficient reduction directly to water (acidic media) or hydroxide ions (alkaline media), whereas values near two indicate predominant hydrogen peroxide formation.

10.8 Hydrogen Peroxide Yield

The percentage of hydrogen peroxide generated during oxygen reduction is calculated as

$$\%H_2O_2 = 200 \cdot (I_R/N) / (I_D + I_R/N)$$

The software calculates and displays the hydrogen peroxide yield as a function of potential and allows comparison of different catalyst materials.

Typical interpretation:

H₂O₂ yield	Interpretation
<5 %	Excellent four-electron catalyst
5–20 %	Mixed reaction pathway
>20 %	Significant peroxide production

10.9 Equilibrium Potential and Reference-Electrode Conversion

For kinetic evaluation, disk potentials are usually expressed as overpotentials relative to the thermodynamic equilibrium (reversible) potential of the reaction under study. SpectroElectroChem Suite can calculate this equilibrium potential automatically from the selected reaction, the pH and the temperature, and convert it to the reference electrode actually used in the experiment.

The pH dependence follows the Nernstian slope $2.303 \cdot R \cdot T / F$ per pH unit (approximately 59 mV per pH unit at 25 °C). The user selects the reaction (for example oxygen reduction, oxygen evolution, hydrogen evolution or hydrogen oxidation), the reference electrode, the pH and the temperature; alternatively a user-defined equilibrium potential can be entered directly.

Scientific check: Verify that the chosen reaction, pH and reference electrode correspond to the real experimental conditions. An incorrect equilibrium potential shifts the entire overpotential axis and therefore the apparent exchange current and Tafel intercept, even though the Tafel slope itself is unaffected.

10.10 Levich and Koutecký–Levich Analysis

The dependence of the diffusion-limited disk current on the rotation rate is evaluated, with the Levich and Koutecký–Levich relations. In the Levich representation the limiting current is plotted against the square root of the angular velocity, $\sqrt{\omega}$. A straight line passing through the origin indicates purely mass-transport-controlled behaviour. In the Koutecký–Levich representation the reciprocal current is plotted against $1/\sqrt{\omega}$; the intercept yields the kinetic current, which is free of mass-transport limitation.

The software determines the limiting current from a user-selected diffusion-limited plateau and reports the plateau drift, the relative scatter and the plateau slope as quality indicators. Absolute cathodic disk currents are used so that reduction currents are handled consistently. Always inspect the selected plateau region before accepting the analysis: a sloping or drifting plateau indicates that a true limiting current has not been reached.

10.11 Tafel Analysis, Exchange Current and Butler–Volmer Parameters

The Tafel module fits the linear region of the potential–log|current| relationship for each rotation rate and reports the Tafel slope, the intercept and the coefficient of determination R^2 . From the intercept and the equilibrium potential the exchange current, or the exchange current density, is derived.

Current density option. Currents may be evaluated either as measured currents or as current densities. When the current density mode is selected, the disk currents are divided by the geometric electrode area (default 0.196 cm², corresponding to a 5 mm diameter disk). Enter the area that matches the electrode actually used.

Mass-transport correction. Optionally, the Tafel region can be corrected for mass transport using the Koutecký–Levich relation, so that the extracted kinetic current more closely reflects the charge-transfer kinetics rather than the combined kinetic and diffusion response. This requires a valid limiting-current plateau (see Section 10.10).

Standard heterogeneous rate constant. If requested, the standard rate constant k_0 is estimated from the exchange current density through the Butler–Volmer relation

$j_0 = n \cdot F \cdot k_0 \cdot c_{\text{O}}^{(1-\alpha)} \cdot c_{\text{R}}^{\alpha}$, using the entered electron number n , transfer coefficient α , and bulk concentrations of the oxidized and reduced species. This calculation is only meaningful for a simple, well-defined one-step redox couple; it is not applied to multi-electron reactions such as ORR, OER, HER or HOR, for which the Butler–Volmer concentration terms are not defined in this simple form.

The Tafel Excel report contains a structured summary sheet with the general settings, the per-curve Tafel results, the limiting-current table (rotation rate, limiting current, drift, relative scatter and plateau slope) and a dedicated section collecting all quality warnings.

10.12 Graphical Output

SpectroElectroChem Suite automatically generates:

- Disk current vs. potential
- Ring current vs. potential
- Combined disk/ring plots
- Hydrogen peroxide yield
- Electron-transfer number
- Interactive three-dimensional RRDE plots
- Excel export of original, background, corrected and processed data

All figures can be exported in publication-quality figure resolution.

Interactive HTML figures allow free rotation, zooming and panning.

By convention, disk currents are reported in milliamperes (mA) and ring currents in microamperes (μA), reflecting the much smaller magnitude of the ring current. This convention is applied consistently in the plots, the combined data sheet and the Excel charts.

Static PNG and PDF graphics are intended for publication-quality figures.

10.13 Best Practice

For reliable RRDE measurements the following recommendations are suggested:

- Polish the disk electrode before each experiment.
- Ensure that the ring electrode is clean.
- Use freshly prepared electrolyte.
- Remove dissolved oxygen when recording background measurements.
- Use identical rotation speeds for all comparative measurements.
- Verify the collection efficiency before quantitative analysis.
- Inspect the limiting current plateau before performing Levich or Koutecký–Levich analysis.

10.14 Worked example: oxygen reduction (ORR) in alkaline solution

This example illustrates a complete RRDE workflow for the oxygen reduction reaction and the values a user would typically enter.

Experimental setup. Glassy-carbon disk carrying the catalyst film and a platinum ring; collection efficiency $N = 0.25$. Electrolyte: 0.1 M KOH saturated with O_2 . The ring is held at a constant potential at which peroxide (HO_2^-) is re-oxidised, for example +0.40 V vs. RHE. The disk is swept from +1.0 V to +0.1 V vs. RHE at 10 mV s⁻¹. Rotation rates: 400, 900, 1600 and 2500 rpm.

Procedure in the software. (1) Import the NOVA CSV file. (2) Subtract an N_2 background for the ring current. (3) Enter the collection efficiency $N = 0.25$. (4) Select the potential window used for the peroxide evaluation. (5) Run the Levich and Koutecký–Levich analysis over the diffusion-limited plateau. (6) Optionally run the Tafel analysis using the ORR equilibrium potential for the chosen pH.

Typical results and interpretation. For a well-behaved four-electron catalyst the disk limiting current scales linearly with $\sqrt{\omega}$ in the Levich plot, the electron number $n(E)$ computed from the disk and ring currents lies close to 4 across the plateau, and the hydrogen-peroxide yield stays below about 5 %. Values of n approaching 2 with a rising % H_2O_2 indicate a two-electron pathway. The Koutecký–Levich intercept gives the kinetic current, and the Tafel slope characterises the charge-transfer kinetics.

11. Signal-processing options

Signal processing is intended to improve data quality while preserving chemically meaningful spectral information. Excessive smoothing or baseline correction may distort quantitative results.

11.1 No processing

Use unprocessed data when the signal quality is sufficient or when the purpose is to document the original measurement. Raw outputs should be retained even when processing is applied.

11.2 Moving average

A moving average replaces each point by a local mean. It is simple and robust but can broaden peaks and shift apparent edges. It is best used for exploratory visualization rather than for quantitative peak-shape analysis.

11.3 Savitzky–Golay filter

The filter preserves polynomial features better than a simple average, but its behaviour depends strongly on the relationship between window length, polynomial order, sampling interval and peak width.

11.4 Baseline subtraction

Baseline methods estimate a low-frequency component. Edge behaviour and broad features require particular care. Record all baseline parameters in reports and retain the calculated baseline in the exported workbook.

11.5 Intensity clipping

Percentile clipping limits extreme values to improve color-scale usage. It should be reported explicitly when figures are used in publications, because it suppresses the displayed magnitude of outliers.

11.6 Logarithmic scaling

log₁₀ scaling increases the visibility of weak features while compressing strong peaks. If negative values are shifted before the transformation, the transformed scale is no longer directly comparable to the original intensity scale.

12. Graphical outputs

12.1 3D surface

The surface plot maps potential, spectral coordinate and response value onto three axes. Rotate and zoom to inspect trends, but confirm interpretations with 2D views, because perspective can hide or exaggerate features.

12.2 Heatmap

The heatmap is often the clearest representation of potential-dependent band evolution. Inspect the color-bar range and check whether clipping or logarithmic scaling has been applied.

12.3 Contour plot

Contours emphasize lines of constant response and can reveal band shifts or transitions. Too many contour levels can suggest a precision that is not present in the measurement.

12.4 Waterfall plot

Waterfall plots show individual spectral traces. They are useful for following peak growth, decay and shift, but the applied offset must be stated or documented.

12.5 Interactive HTML

HTML files open in a web browser. The Plotly controls allow rotation, zooming, panning and image export. Some files may load the Plotly JavaScript library from a content-delivery network; offline use may therefore depend on how the HTML file was generated.

13. Excel and file exports

The exported Excel workbooks are intended to provide full numerical traceability of all graphical results.

13.1 Output location

If an output folder is selected, all generated files are written there. Otherwise, many modules use the folder containing the source CSV file. Use a dedicated analysis folder to avoid mixing raw and processed data.

13.2 Metadata

Metadata sheets document the source file, the selected range, the matrix dimensions, the smoothing method and the scaling method. Keep this sheet with the numerical outputs when data are transferred to collaborators.

13.3 Matrix and long formats

- Matrix format is convenient for surface and heatmap creation.
- Long format contains one row per spectral-coordinate/potential pair and is convenient for statistical software and database import.

13.4 File naming

Generated file names commonly include the source stem, the selected spectral range, the smoothing method, the scaling method and the plot type. Do not rename files in a way that removes processing information unless the information is preserved elsewhere.

14. Recommended scientific workflow

- Archive the original instrument export as read-only raw data.
- Create a working copy in CSV format.

- Check the axis orientation, units, missing values and matrix dimensions.
- Run the analysis first with no smoothing and no clipping.
- Inspect the raw plots for spikes, discontinuities and import errors.
- Apply the least aggressive processing that adequately addresses the analytical problem.
- Inspect the baseline and smoothing control plots.
- Export all numerical stages, not only the final figure.
- Record the software version, DOI and all parameters in the laboratory notebook.
- For publication, state whether values were baseline-corrected, smoothed, clipped, logarithmically transformed or vertically offset.
- Archive both raw data and processed results.

15. Troubleshooting

Problem	Likely cause	Recommended action
Batch file opens and closes immediately	Python or a package is missing	Open Command Prompt in the program folder and run the batch file there to read the error message.
“python” is not recognized	Python is not on PATH	Reinstall Python with “Add Python to PATH” enabled, or edit the batch file to use the py launcher.
ModuleNotFoundError	Packages were installed into another Python environment	Run “python -m pip install -r requirements.txt” using the same interpreter that starts the program.
CSV file not accepted	Non-numeric cells, wrong delimiter, incorrect matrix orientation or encoding	Open the file in a text editor, verify the first row/column, and export it again as CSV UTF-8.
No data points in the selected range	Range lies outside the detected axis values	Use the range-check function and enter limits within the available interval.
Savitzky–Golay error	Window is even, too	Use an odd window larger than the polynomial order and smaller than the number of spectral points.

Problem	Likely cause	Recommended action
	small or longer than the available data	
Plot dominated by a single spike	Cosmic ray or extreme outlier	Inspect the raw data; correct the raw measurement if scientifically justified, or use percentile clipping for visualization only.
Baseline removes real bands	Baseline parameters are too aggressive	Reduce the polynomial flexibility or penalty settings and inspect the baseline overlay.
Interactive HTML is blank offline	Plotly JavaScript is referenced online	Connect to the Internet, or regenerate the HTML with embedded Plotly JavaScript.
Excel file cannot be overwritten	The workbook is open in Excel	Close the workbook, or select a different output file name or folder.
Application does not start from another folder	An old batch file used a fixed path	Use the release batch file, which starts relative to its own location, and keep the project files together.
Memory Error during export	Very large matrices	Reduce the spectral range or process the dataset in smaller sections

If interactive HTML plots do not open automatically, verify that a default web browser is configured in Windows.

16. Reproducibility and data integrity

- Retain the original CSV file and the instrument-native file.
- Record the exact software version and the version-specific DOI.
- Record all processing parameters and the selected spectral ranges.
- Record all processing parameters together with the software version and DOI.
- Preserve the exported metadata, raw matrices and processed matrices.
- Do not infer quantitative concentrations from arbitrary intensity values without calibration.
- Use the same processing settings when comparing samples, unless a justified exception is documented.

- Verify that the exported figures correspond to the numerical workbook used in the publication.
- Archive the generated Excel workbook together with the corresponding HTML and PDF outputs.

Version-specific citation: For analyses performed with version 5.5.3, cite DOI 10.5281/zenodo.21283231. The concept DOI 10.5281/zenodo.21283231 refers to the software project across versions and resolves to the latest archived release.

17. Citation and licensing

17.1 Recommended citation

Users are encouraged to cite both the software DOI and any associated scientific publications where appropriate.

Habekost, A. (2026). SpectroElectroChem Suite (Version 5.5.3) [Computer software]. Zenodo. <https://doi.org/10.5281/zenodo.21283231>

17.2 Software availability statement

SpectroElectroChem Suite is open-source software for the analysis and visualization of Raman, SERS, absorption and fluorescence spectro-electrochemical data. The source code is available from the GitHub repository, and the version used in this work is permanently archived on Zenodo (DOI: 10.5281/zenodo.21283231).

17.3 License

The software is distributed under the MIT License. Users may use, copy, modify and redistribute the software subject to the license terms. The software is provided without warranty. Consult the LICENSE file included in the repository for the full legal text.

17.4 Acknowledgement

Parts of the source code were developed with the assistance of OpenAI ChatGPT and were subsequently reviewed, scientifically validated, modified and substantially extended by the author.

Appendix A. CSV templates

A.1 Raman/SERS voltammogram template

	-0.40	-0.30	-0.20	-0.10
300	120	125	131	140
305	118	124	130	139
310	119	126	134	143
315	123	129	138	148

A.2 Absorption/fluorescence voltammogram template

	0.00	0.10	0.20	0.30
400	0.021	0.028	0.039	0.051
405	0.023	0.031	0.043	0.056
410	0.026	0.035	0.048	0.062
415	0.030	0.040	0.054	0.069

A.3 Raman spectrum template

300	120
305	118
310	119
315	123

A.4 Absorbance spectrum template

300	0.120
305	0.118
310	0.119
315	0.123

A.5 Fluorescence spectrum template

300	120
-----	-----

305	118
310	119
315	123

Appendix B. Output-file glossary

Output	Meaning
Raw_Intensity_Matrix	Imported numerical matrix after axis and empty-cell validation.
Smoothed_Intensity_Matrix	Matrix after the selected smoothing operation.
Plotted_Intensity_Matrix	Values actually used for the visual output, after optional clipping or logarithmic scaling.
Waterfall_Unshifted_Values	Processed spectra before vertical displacement.
Waterfall_Shifted_Values	Spectra after addition of the waterfall offsets.
Waterfall_Offsets	Numerical offset applied to each spectrum.
Metadata	Source, range, dimensions and processing settings.
*.html	Interactive browser-based plot.
*.png	Static raster figure.
*.pdf	Static vector or page figure, depending on the export method.
*.xlsx	Excel workbook with numerical data and metadata.

Appendix C. Quick-start checklist

- Download the official ZIP archive from Zenodo or GitHub.
- Extract the ZIP archive completely.
- Install Python 3.10 or newer.
- Run Install_required_Python_packages.bat once.
- Run Test_Installation.bat.
- Start the application with Start_SpectroElectroChem_Suite.bat.
- Alternatively, start the ready-to-run Windows executable (SpectroElectroChem_Suite.exe) or use the installed shortcut – no Python installation is required in that case.
- Select the correct analysis module.
- Load a correctly formatted CSV file.

- Choose a valid spectral range.
- Use conservative processing parameters.
- Run the analysis and inspect the control plots.
- Keep the raw, processed and metadata outputs together.
- Cite version 5.5.3 with DOI 10.5281/zenodo.21283231.

Appendix D. Typical RRDE Parameters

Parameter	Typical value
Collection efficiency	0.25
Rotation	500–3000 rpm
O ₂ concentration	1 mmol/L
Potential range	0.7 V to -1.5 V

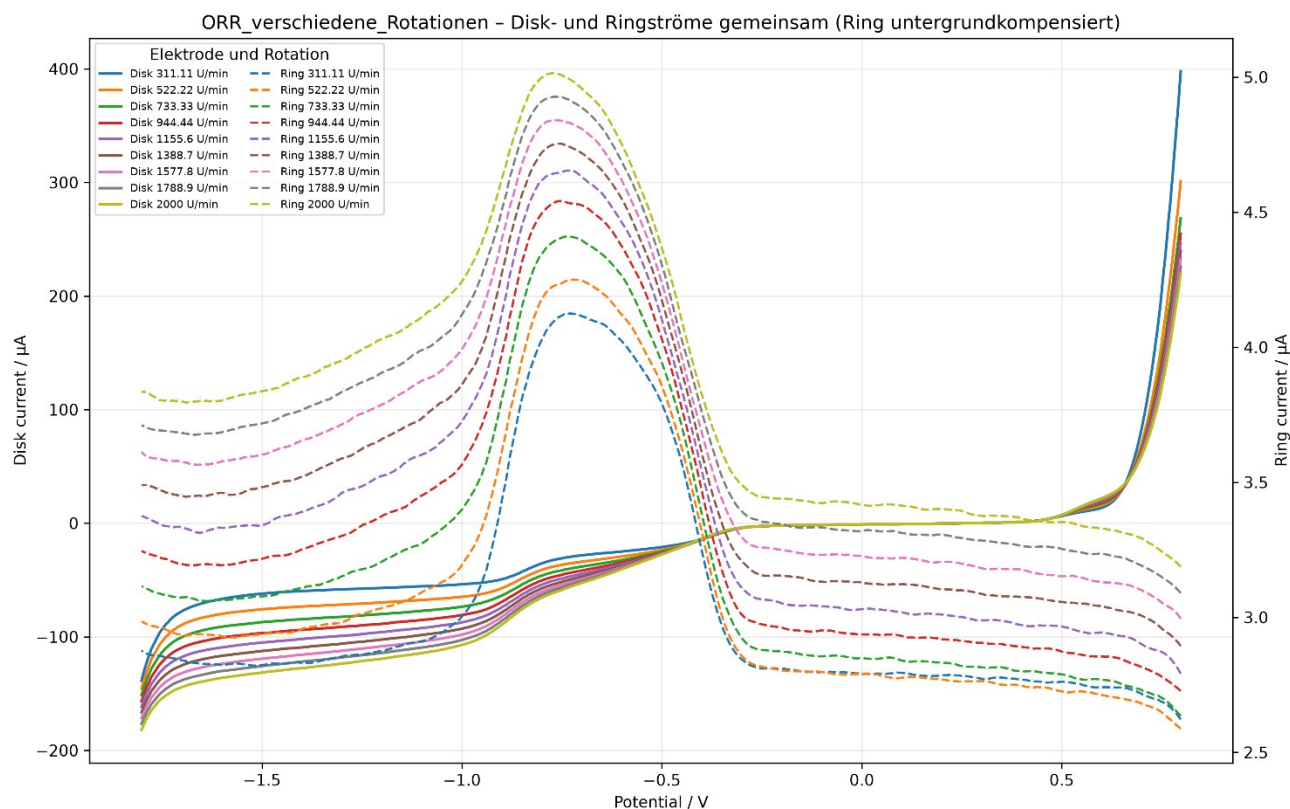


Figure 8: RRDE results for the oxygen reduction reaction (ORR): disk and ring currents between 0.7 and -1.8 V.

Appendix E. Installation from Python source (alternative route)

This appendix condenses the Python-source installation for advanced users and developers. Most users should instead use the executable or installer described in Section 4.5.

- 1. Install Python.** Install Python 3.10 or newer and enable “Add Python to PATH” during setup.
- 2. Install packages.** In the extracted program folder, run `Install_required_Python_packages.bat` (or `py -m pip install -r requirements.txt`).
- 3. Test.** Run `Test_Installation.bat` to verify that Python and the principal dependencies import correctly.
- 4. Start.** Launch the suite with `Start_SpectroElectroChem_Suite.bat`, or from a Command Prompt with `python main.py`.